


# TheraCell



TheraFuze DBF®

FIBER BULLETS™

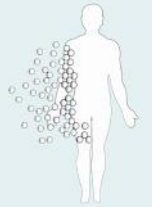

FIBER  
BULLETS™

6mm CALIBER FIBER BULLETS IN PACKS OF 5 OR 10

TheraFuze DBF®

Demineralized Bone Fiber Bullet Implants for MIS Delivery created with Osteoconductive Demineralized Bone Fibers with Osteoinductive Potential.

5 Pack of 6 mm Caliber BULLETS

ORDER: TC-FBL-0605



10 Pack of 6 mm Caliber BULLETS

ORDER: TC-FBL-0610



<b>What is it?</b>	A demineralized bone fiber (DBF) form designed for use in minimal access surgery.
<b>What is it made of?</b>	Demineralized bone fibers molded into a proprietary shape that can be delivered through a 6mm i.d. or larger tube. 100% cortical bone, with no excipients. Produced in compliance with HCT/P regulations under Title 21 Code of Federal Regulations (CFR), Part 1271.
<b>How can it be used?</b>	The Fiber Bullets™ can be inserted into a bone void or cage through a 4mm or a 6mm i.d. tube delivery system and will hydrate in the site. The lyophilized Fiber Bullets™ can be delivered using minimally invasive techniques into the desired site delivering osteoconductive demineralized bone fibers with osteoinductive potential to stimulate new bone formation for long term structural stability.
<b>How does it work?</b>	The Fiber Bullets™ swell as they hydrate, forming an entangled mass that resists migration and enhances the ability to fill defects completely, e.g. to achieve endplate to endplate contact, which is important to achieving fusion.
<b>What sizes are available?</b>	6mm Caliber Fiber Bullets in packs of 5 or 10.
<b>Patents</b>	The implants and/or techniques associated with this product may be covered by one or more of the United States Patents: US 9,572,912 US 9,486,557 AUS 2014253753 Other patents pending

TO ORDER OR FOR MORE INFORMATION PLEASE CONTACT US AT - 630-953-9594

**THERACELL, INC.**

14930 VENTURA BLVD., SUITE #325  
SHERMAN OAKS, CA 91403

630-953-9594

[www.theracellinc.com](http://www.theracellinc.com)

TheraFuze DBF® is a registered trademark of TheraCell, Inc.  
Fiber Bullets™ is a trademark of TheraCell, Inc.

ALL RIGHTS RESERVED.

