


TheraCell



TheraFuze DBF®

FIBERFORM™+

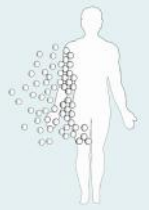

FIBERFORM™+

1cc, 5cc & 10cc MOLDABLE DBF WITH MINERALIZED CORTICAL CHIPS

TheraFuze DBF®

A moldable demineralized bone fiber form designed for use in orthopedic and spine surgery that also includes mineralized cortical chips.



10 cc
ORDER: TC-FFP-10



5 cc
ORDER: TC-FFP-05



1 cc
ORDER: TC-FFP-01

What is it?	A moldable demineralized bone fiber (DBF) form designed for use in orthopedic and spine surgery that also includes mineralized cortical chips.
What is it made of?	Demineralized bone fibers plus added mineralized cortical bone chips. 100% cortical bone, with no excipients. Produced in compliance with HCT/P regulations under Title 21 Code of Federal Regulations (CFR), Part 1271.
How can it be used?	Can be rehydrated with saline, blood or bone marrow aspirate.
How does it work?	The hydrated FiberForm™ + is a conformable excipient-free putty that can be placed in any bone defect or area where new bone formation is required. Cortical chips provide additional compression resistance and radiodensity for visualization. Osteoconductive with osteoinductive potential to stimulate long term bone formation.
What sizes are available?	1 cc 5 cc 10 cc
Patents	The implants and/or techniques associated with this product may be covered by one or more of the United States Patents: US 9,572,912 US 9,486,557 AUS 2014253753 Other patents pending

TO ORDER OR FOR MORE INFORMATION PLEASE CONTACT US AT - 630-953-9594

THERACELL, INC.
14930 VENTURA BLVD., SUITE #325
SHERMAN OAKS, CA 91403

630-953-9594

www.theracellinc.com

TheraFuze DBF® is a registered trademark of TheraCell, Inc.
FiberForm™ is a trademark of TheraCell, Inc.

ALL RIGHTS RESERVED.

