


TheraCell



TheraFuze DBF[®]

FIBER WRAP[™]

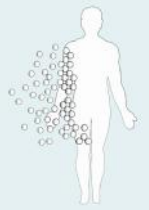

FIBER WRAP[™]

5 cm x 5.5 cm x 1mm • 10 cm x 5.5 cm x 1 mm • 10 cm x 11 cm x 1mm

TheraFuze DBF[®]

Demineralized Bone Fiber Sheet with FiberLok™ in a Versatile Graft Containment Form



10 cm x 11 cm x 1mm
ORDER: TC-FW-101101



10 cm x 5.5 cm x 1mm
ORDER: TC-FW-105501



5 cm x 5.5 cm x 1mm
ORDER: TC-FW-055501

What is it?	Demineralized bone fibers molded into a 1mm thick sheet with FiberLok™ technology.
What is it made of?	100% cortical bone, with no excipients. Produced in compliance with HCT/P regulations under Title 21 Code of Federal Regulations (CFR), Part 1271.
How can it be used?	The lyophilized sheet may be rehydrated with saline or the patient's blood or bone marrow creating a flexible sheet that can be manipulated at the surgeon's discretion. For example: <ul style="list-style-type: none"> • As a versatile lower profile alternative to boats for PLF procedures. • Rolled up with autograft to make a "cigar" that can be laid in the posterolateral gutters for PLF. • Periosteum replacement in trauma applications to contain bone fragments / graft. • TJR (acetabular reconstruction). • As a bone bandage.
How does it work?	TheraCell's long and strong DBF fibers provide a contiguous matrix that is enhanced with our proprietary FiberLok™ technology that provides strength to the hydrated sheet to allow it to be manipulated after rehydration by the surgeon. Osteoconductive with osteoinductive potential to stimulate long term bone formation.
What sizes are available?	10 cm x 11 cm x 1mm ● 10 cm x 5.5 cm x 1 mm ● 5 cm x 5.5 cm x 1mm
Patents	The implants and/or techniques associated with this product may be covered by one or more of the United States Patents: US 9,572,912 US 9,486,557 AUS 2014253753 Other patents pending

TO ORDER OR FOR MORE INFORMATION PLEASE CONTACT US AT - 630-953-9594

THERACELL, INC.

14930 VENTURA BLVD., SUITE #325
SHERMAN OAKS, CA 91403

630-953-9594

www.theracellinc.com

TheraFuze DBF® is a registered trademark of TheraCell, Inc.
Fiber Wrap™ and FiberLok™ are trademarks of TheraCell, Inc.

ALL RIGHTS RESERVED.

